




STARMed CaStar R

Hood for non-invasive mechanical ventilation (NIV)

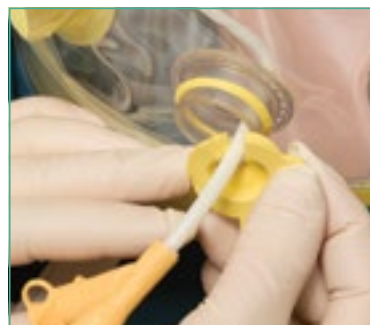

Critical Care ▪ Patient Interfaces ▪ **STARMed**



STARMed CaStar R hood for NIV

StarMed CaStar R

The CaStar R is a comfortable, versatile and lightweight patient interface for non-invasive mechanical ventilation (NIV) designed to offer an alternative to traditional methods of NIV delivery. The CaStar R is well tolerated for extended periods and can reduce the risks associated with both endotracheal intubation and face mask use. It can be used in semi-recumbent, supine and prone positions providing an ideal option for NIV therapy.

Technical Features

- Manufactured from phthalate-free materials
- 22M connectors compliant to ISO standard
- Comfortable and atraumatic underarm fastening straps
- Built in inflatable neck cushion

- Two versatile sealed access ports suitable for single or multiple lumen probes or catheters
- New design patient access port. Fitted with bi-directional anti-asphyxiation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated

Bibliography:

- Mojoli F, et al. An optimized set-up for helmet noninvasive ventilation improves pressure support delivery and patient-ventilator interaction, *Intensive Care Med.* 2013 Jan;39(1):38-44
- Vargas F, et al. Helmet with specific settings versus facemask for noninvasive ventilation, *Crit Care Med.* 2009 Jun;37(6):1921-8
- Mojoli F, et al. Carbon dioxide rebreathing during non-invasive ventilation delivered by helmet: a bench study, *Intensive Care Med.* 2008 Aug;34(8):1454-60
- Antonelli M, et al. New treatment of acute hypoxemic respiratory failure: Noninvasive pressure support ventilation delivered by helmet - A pilot controlled trial, *Crit Care Med.* 2002 Mar;30(3):602-8



CaStar R hood for NIV with inflatable neck cushion, bi-directional anti-asphyxiation valve and patient access port*

Code	Size	Neck circumference	Box Qty.
CP211XXL/2R	Extra Extra Large	≥50cm	5
CP211XL/2R	Extra Large	45-52cm	5
CP211L/2R	Large	40-47cm	5
CP211M/2R	Medium	34-41cm	5
CP211S/2R	Small	27-34cm	5
CP211XS/2R	Extra Small	17-27cm	5



www.intersurgical.nl/info/starmed

Video available at www.intersurgical.nl

Make an enquiry

*Tubing not included and has to be ordered separately. Please refer to the website for the range of accessories.



Intersurgical Benelux B.V. Vluchttoord 28, 5406 XP, Uden, Nederland

T: +31 (0)413 243 860 F: +31 (0)413 243 869 info@intersurgical.nl www.intersurgical.nl



Manufacturer:
Intersurgical S.p.A.
Mirandola (MO) Italy

Please think before you print

Save energy and paper. If you must print this information sheet please print it double sided.

IS1.18_NL ▪ Issue 8 06.20

UK ▪ Ireland ▪ France ▪ Germany ▪ Spain ▪ Portugal ▪ Italy ▪ Benelux ▪ Sweden ▪ Denmark ▪ Lithuania ▪ Russia ▪ Czech Republic
Turkey ▪ South Africa ▪ China ▪ Japan ▪ Taiwan ▪ Philippines ▪ USA ▪ Canada ▪ Colombia ▪ Australia