



FaceFit™ NIV masks

Lightweight and comfortable NIV (Non-invasive Ventilation) masks for conventional ventilation or CPAP and bilevel applications



Critical Care ▪ Patient Interfaces



FaceFit™ NIV masks

FaceFit is a comfortable single use adult hospital mask for use in non-invasive ventilation.

Two options are available; one option is for use with conventional dual limb breathing systems, and the second for bilevel and CPAP applications where an anti-asphyxiation valve and CO₂ leak port is required.

FaceFit masks are available in small, medium and large sizes to fit a wide range of patients. A size guide is included to help with making the right choice prior to opening the packaging.

Adjustable rear head pad

Allows for optimal fit for a wide range of patients

Lightweight

Improving patient tolerance

Easily adjustable harness and quick release mechanism

Allowing for fast removal of the mask from the patient and optimum fit

Video available at www.intersurgical.nl



Anti-asphyxiation valve (FaceFit ported only)

Prevents rebreathing in the case of inadequate ventilation

CO₂ Leak port (FaceFit ported only)

For effective CPAP or bilevel ventilation

Low dead space


For efficient respiratory support

6mm stem connector with retainable cap

For supplementary oxygen or monitoring

Soft silicone seal

Ensuring optimum performance

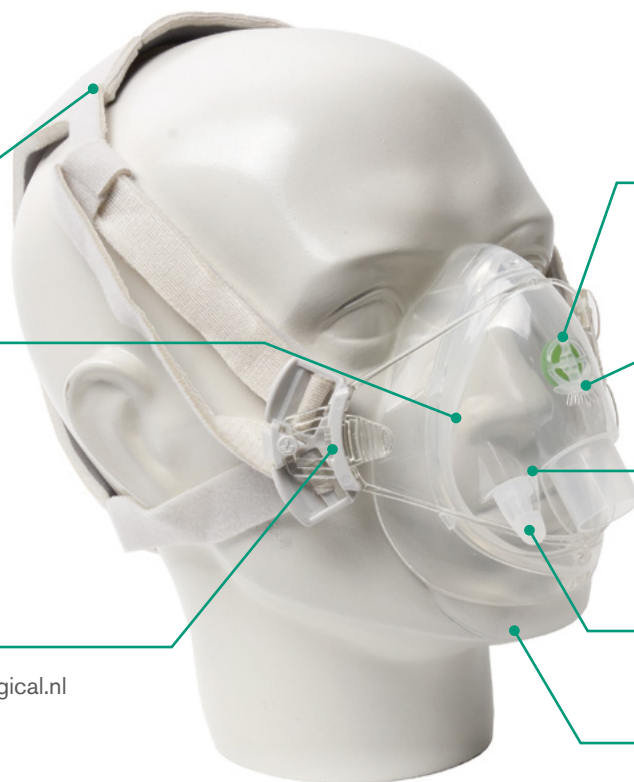

FaceFit NIV Mask - without anti-asphyxiation valve and without CO₂ leak port

For use with dual limb breathing systems

illustrated

Code	Description	Box Qty.
2257000	FaceFit, NIV mask, large adult	12
2256000	FaceFit, NIV mask, medium adult	12
2255000	FaceFit, NIV mask, small adult	12

FaceFit NIV Ported Mask - with anti-asphyxiation valve and CO₂ leak port

For use with single limb breathing systems that do not have a CO₂ leak port

Code	Description	Box Qty.
2252000	FaceFit, NIV Ported Mask, with anti-asphyxiation valve and CO ₂ leak port, large adult	12
2251000	FaceFit, NIV Ported Mask, with anti-asphyxiation valve and CO ₂ leak port, medium adult	12
2250000	FaceFit, NIV Ported Mask, with anti-asphyxiation valve and CO ₂ leak port, small adult	12

Make an enquiry



Intersurgical Benelux B.V. Vluchttoord 28, 5406 XP, Uden, Nederland

T: +31 (0)413 243 860 F: +31 (0)413 243 869 info@intersurgical.nl www.intersurgical.nl



The manufacturer Intersurgical Ltd is certified to ISO 9001:2008, ISO 13485:2003 and ISO 14001:2015

Please think before you print

Save energy and paper.
If you must print this information sheet please print it double sided.

IS1.15 NL • Issue 6 03.18