



Laryngoscope blades


For single patient use combining safety and strength



Airway Management ▪ Airway Devices

Plastic single use laryngoscope blades

Intersurgical plastic single use laryngoscope blades provide a safe, strong, light and cost-effective alternative to reusable metal versions. The range of Miller and Macintosh blades provides effective illumination and visualisation of the airway to aid with intubation and airway clearance. Our blades are compatible with the ISO7376-3 green specification handles, as indicated by the green dot near the heel of the blade.

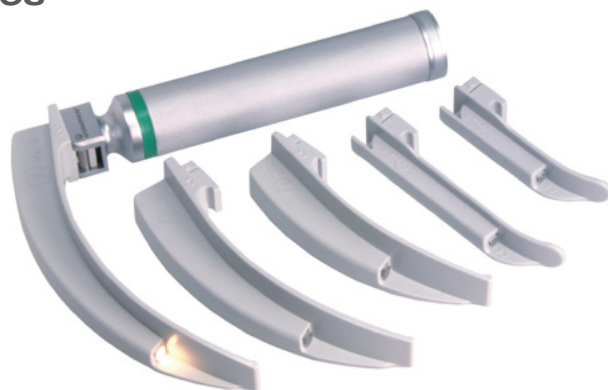

High optical light transmission

Provides the user with an excellent light transfer to the blade tip

Single construction blade made from high impact grade ABS plastic

Rigid blade which is a cost effective, disposable alternative to reusable and disposable metal blades

Engagement to handle click mechanism via two plastic lugs

Provides secure and positive connection of the blade to the handle



Ultra smooth blade with matt finish

Smooth transition into the airway and easy passage of the tube, with no light reflection off the blade

Macintosh style



Code	Description	Size	Box Qty.
7044000	Plastic laryngoscope blade, fibre optic, Macintosh	4	40
7043000	Plastic laryngoscope blade, fibre optic, Macintosh	3	40
7042000	Plastic laryngoscope blade, fibre optic, Macintosh	2	50

Miller style



Code	Description	Size	Box Qty.
7041000	Plastic laryngoscope blade, fibre optic, Miller	1	70
7040000	Plastic laryngoscope blade, fibre optic, Miller	0	70

Make an enquiry

Please think before you print

Save energy and paper. If you must print this information sheet please print it double sided.

IS5.10 ▪ Issue 8 01.18



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
 T: +44 (0)118 9656 300 F: +44 (0)118 9656 356 info@intersurgical.com www.intersurgical.com

UK ▪ Ireland ▪ France ▪ Germany ▪ Spain ▪ Portugal ▪ Italy ▪ Benelux ▪ Sweden ▪ Denmark ▪ Lithuania ▪ Russia ▪ Czech Republic
 Turkey ▪ South Africa ▪ China ▪ Japan ▪ Taiwan ▪ Philippines ▪ USA ▪ Canada ▪ Colombia ▪ Australia