

QuadraLite™

a truly **unique** anaesthetic face mask



QuadraLite™

a truly **unique** anaesthetic face mask

Soft seal

A comfortable airtight fit with flexibility for different face sizes

Transparent shell
For excellent visibility

Colour coding

Colour matching with other Intersurgical products for easy selection

Anatomical shape

Provides a precise fit combined with lower dead space

Unique Design Developed by Intersurgical

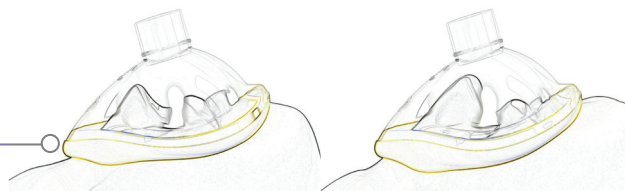
QuadraLite™ does not use an inflatable cushion. Innovative moulding techniques have enabled the combination of two ideally suited materials allowing the creation of a face mask with an anatomical shape and soft seal. This provides a comfortable, superior airtight fit with a crystal clear shell to enhance patient visibility.

Four Sizes for all Patients





The mask is lightweight and requires only minimal pressure to achieve a seal over a wide range of face sizes and shapes. Only four sizes instead of the traditional seven provide a solution for all patients.

The Greener, Safer Choice

QuadraLite™ is part of the Intersurgical Eco Range of products containing no PVC which are designed to have a low environmental impact. Masks are phthalate free to protect the patient and a serial number for patient records and traceability is etched directly onto every mask.

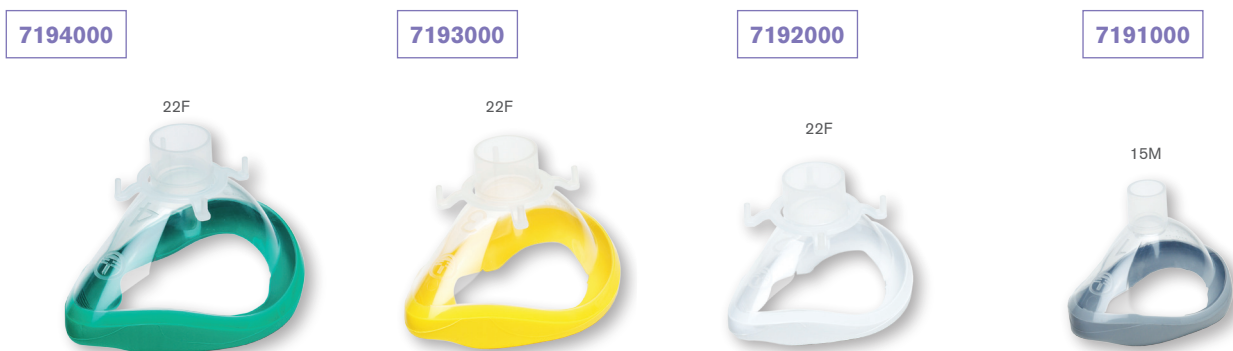







Our quick sizing guide for comparison

Patient	Age	Weight (Kg)	Traditional mask size	QuadraLite™ masks
Neonatal	Preterm	< 2	0	 <p>Code: 7191000</p>
	Term	2 – 4		
Infant	3 months	6	1	
	6 months	8		
	9 months	9		
Paediatric	1 year	10	2	
	2 years	12		
	3 years	14		
	4 years	16		
	5 years	18		
Small Adult	6 years	20	3	 <p>Code: 7193000</p>
	7 years	23		
	8 years	26		
	9 years	28		
	10 years	32		
	11 years	35		
Medium Adult	12 years	39	4	
	13 years	44		
	14 years	50		
	15 years	54		
Large Adult	16 years	58	5	 <p>Code: 7194000</p>
	16+ years	60 – 120		
Extra Large Adult	16 + years	> 120	6	 <p>Code: 7194000</p>

Ordering information

QuadraLite™ anaesthesia mask



Code	7194000	7193000	7192000	7191000
Box qty.	35	35	25	30
Mask size	5-6	3-4	2-3	0-1
Patient size	Large	Medium	Small	Infant
Connection size	22F	22F	22F	15M
Colour code				
Hook ring	✓	✓	✓	
	✓	✓	✓	✓

Code	Product description	Box qty.
7089000	Anaesthetic mask hook ring for both QuadraLite™ and EcoMask™ with 22mm port	30

ISS.24 NL • Issue 2 10.13