


Dual Lumen Nasal Cannulae

For use with Portable Liquid Oxygen Systems which conserve oxygen by supplying gas only when the patient demands

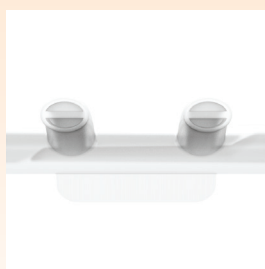

Dual Lumen Nasal Cannulae

Features and Benefits

Multi channel kink resistant tubing

Allows oxygen to continue to flow if the tubing becomes kinked

Satin material
20% less friction¹

Interconnected tubing
Effective tube management

Delivery & sensing channels in both cannulae prongs

optimal demand readings and avoids imbalance between nostrils


Soft prongs
Minimise irritation

The Dual Lumen cannulae is specifically designed for use with all electronic or pressure activated portable liquid oxygen systems such as the HELiOS™. It features a dual lumen cannulae which delivers oxygen to, and senses demand from both nostrils. When demand drops so does the supply of oxygen thus conserving the available supply of gas. The cannulae is made from Intersurgical's soft satin material, proven in testing to create 20% less friction¹ to the patients face and ears compared to a standard cannulae.

Delivers Oxygen
Senses pressure

Ordering information

1165011

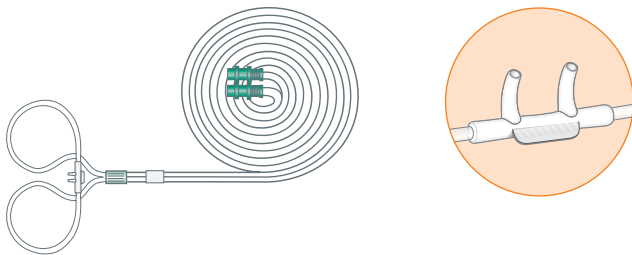

| | | | |
|---------|--|----|--|
| 1165011 | Dual Lumen, adult, nasal cannulae with curved prongs and tube, 1.25m | 50 | |
| 1165012 | Dual Lumen, adult, nasal cannulae with curved prongs and tube, 1.8m | 50 | |
| 1165013 | Dual Lumen, adult, nasal cannulae with curved prongs and tube, 2.5m | 50 | |

HELiOS is a registered trademark of CAIRE® Inc.

References

1. Annan. J. 69381 15.04.13

IS10.30 NL • Issue 1 10.13



Intersurgical Ltd,
Crane House, Molly Millars Lane,
Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 9656 300 F: +44 (0)118 9656 356
info@intersurgical.com www.intersurgical.com

Deutschland
info@intersurgical.de

France
info@intersurgical.fr

Italia
info@intersurgical.it

España
info@intersurgical-es.com

Portugal
info@intersurgical.pt

Nederland
info@intersurgical.nl

Lietuva
info@intersurgical.lt

Россия
info@intersurgical.ru

Česká Republika
info@intersurgical.cz

Taiwan
info@intersurgical.com.tw

South Africa
info@intersurgical.co.za

Philippines
info@intersurgical.ph

Japan
info@intersurgical.co.jp

USA
support@intersurgicalinc.com

China
info@intersurgical-cn.com

