



# Intersurgical EcoLite™ adult high concentration oxygen mask



Oxygen and Aerosol Therapy ▪ Oxygen Therapy



Kwaliteit, innovatie en keuze



@intersurgicalNL



/intersurgical NL

www.intersurgical.be



## Intersurgical EcoLite™ adult high concentration oxygen mask

### Welcome to the Comfort Zone

The unique design and materials used to manufacture the Intersurgical EcoLite mask provides improved patient comfort and a reduced environmental impact.

### Comfortable for the Patient

Improved patient comfort has been key to the development of the Intersurgical EcoLite mask. The latest manufacturing technology has enabled us to combine two non-PVC materials in the same product, the polypropylene material forming the body of the Intersurgical EcoLite mask is clear, lightweight and rigid enough to maintain the masks shape, whilst the second softer TPE material is utilised in the seal, which is in contact with the patients face.

### Comfortable for the Environment

The Intersurgical EcoLite mask is an important part of our Eco range, designed as part of our ongoing focus on sustainable development. The use of PVC in medical products has been questioned and its impact on the environment. The utilisation of alternative materials has resulted in the elimination of PVC from the Intersurgical EcoLite mask, reducing the environmental impact.

### Intersurgical EcoLite

Further information available at [www.intersurgical.be/info/IntersurgicalEcoLite](http://www.intersurgical.be/info/IntersurgicalEcoLite)

### Eco range

Further information available at [www.intersurgical.be/info/Eco](http://www.intersurgical.be/info/Eco)



### Incurved nose seal

Conforms to different nose shapes designed to minimise the risk of oxygen entering patient's eyes

### No metal nose clip

MRI compatible

### Clear mask

permits visualisation of the patient's nose and lips

### Two integrated safety valves

Allows patient to breathe room air in the event oxygen supply failure

### Soft face seals



Increased patient comfort

### Elastic can be positioned under or over the ears

Below ear position eliminates trauma to top of ears

### A choice of "under chin" or "on chin" positions

Provides a better fit on a wider range of patient face shapes

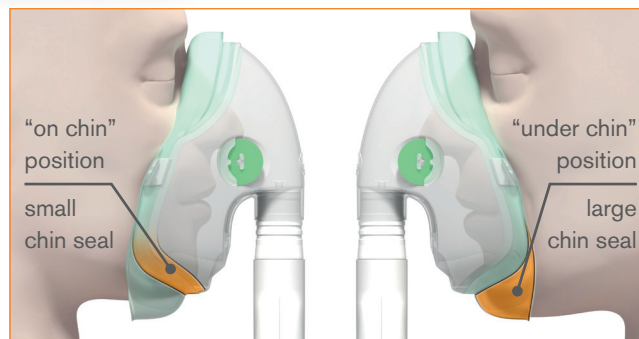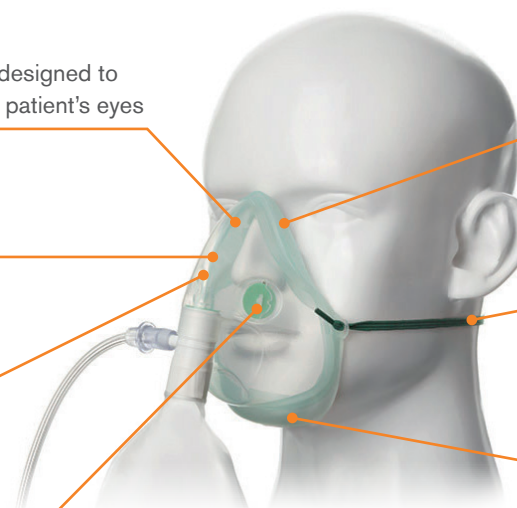

Code	Description	Tube length	Box Qty.
1181015	Intersurgical EcoLite, adult, high concentration oxygen mask with tube	2.1m	24

Video available at [www.intersurgical.be/info/videos](http://www.intersurgical.be/info/videos)

Lower environmental impact product. Note: reservoir bag and tube contain PVC.

IS10.17 BENL • Issue 9 12.22



Intersurgical Benelux B.V. Vluchttoord 28, 5406 XP, Uden, Nederland  
T: +31 (0)413 243 860 F: +31 (0)413 243 869 [info@intersurgical.be](mailto:info@intersurgical.be) [www.intersurgical.be](http://www.intersurgical.be)



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016, ISO 14001:2015 and MDSAP



UK • Ireland • France • Germany • Spain • Portugal • Italy • Benelux • Sweden • Denmark • Lithuania • Russia • Czech Republic  
Turkey • South Africa • China • Japan • Taiwan • Philippines • USA • Canada • Colombia • Australia